


# ***Orthanc-Tools-Js***

*Documentation Guide*

***Salim Kanoun***

*Salim.kanoun@gmail.com*

# I. Architecture

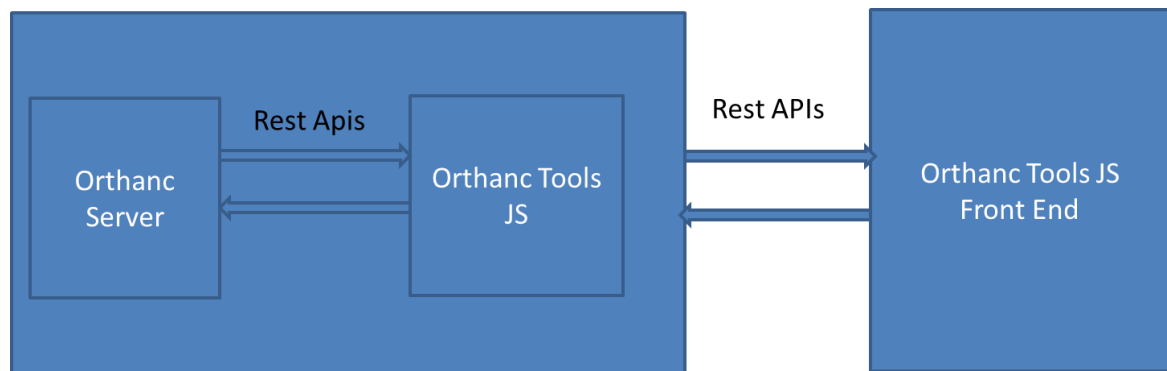

Orthanc Tools JS software is built on the top of Orthanc and extending its capabilities by:

- Providing a super user-friendly interface to perform complex DICOM operations
- Adding feature to Orthanc such as
  - o Role based access control
  - o Automation in Retrieve, Anonymization to perform batches of DICOM operations
  - o Integration of third party app such as viewers like OHIF or StoneOfOrthanc

Technically, OrthancToolsJS create a Restful API server on the top of Orthanc, this server extends Orthanc capabilities of Orthanc by adding additional high level APIs for some features while it act as a simple reverse proxy for other DICOM features that fully rely on Orthanc built'in features.

The built-in web Front end which is a React web app that can be accessed with any browser, the backend is written on NodeJS/ExpressJS.

OrthancToolsJS licence is Affero General Public Licence v3. This license guarantees the free use of this software for any purpose, without any restrictions (neither time nor number of installations).

Pay attention that according to the license terms

- No warranty are associated with the software
- Modification of the source code **must** be redistributed under the same license and this even if your changes are for your own usage only.
- Software solutions that will be dependent of OrthancToolsJS are prone to be affected by this copyleft licence and **must** also disclose its source code.

More information about licensing is available at <https://www.fsf.org/>

At this stage we do not provide offers of services for OrthancToolsJS, however if you have needs of support service or need of a license exception feel free to contact the maintainer ([salim.kanoun@gmail.com](mailto:salim.kanoun@gmail.com))

## II. Installation & upgrade

### a. With Docker

Orthanc deployment is preferred as it ships a full install ready to use (NodeJS + OrthancToolsJs)

The Docker Image name is in DockerHub: salimkanoun/orthanc-toolsjs

Orthanc Tools JS needs the following software to run

- HTTPS reverse proxy
  - o Add the X-Forwarded-For header for client ip logs
- Redis
- A database (soon : Postgres or Mysql)
- Orthanc (obviously)

A docker compose sample for deployment of Orthanc and OrthancToolsJs is available here: <https://github.com/salimkanoun/Orthanc-Tools-JS/blob/master/docker-compose.yml>

See the docker-compose to set the environment variable

To Benefit for AET/Peers declaration from OrthancToolsJS and to benefit for OHIF integration, edit the following settings in orthanc JSON

- ORTHANC\_\_DICOM\_MODALITIES\_IN\_DATABASE: "true"
- ORTHANC\_\_ORTHANC\_PEERS\_IN\_DATABASE: "true"
- *ORTHANC\_\_DICOM\_WEB\_\_ENABLE: "true"*
- *ORTHANC\_\_DICOM\_WEB\_\_ROOT: "/dicom-web/"*
- *ORTHANC\_\_DICOM\_WEB\_\_ENABLEWADO: "true"*
- *ORTHANC\_\_DICOM\_WEB\_\_WADOROOT: "/wado"*
- *ORTHANC\_\_DICOM\_WEB\_\_STUDIESMETADATA: "MainDicomTags"*
- *ORTHANC\_\_DICOM\_WEB\_\_SERIESMETADATA: "Full"*

*Italicized items are for OHIF integration*

### b. Upgrade

Until we will reach a stable version >1.0.0 the database schema will not be stable nor upgradable, so you will need to wipe all previous install and do a fresh install

### c. Other OS

It would be theoretically possible to run Orthanc Tools JS in Windows, Linux, MacOS without docker.

However we choose to not support these installs as it will rely on unofficial ports (such as windows version of Redis).

If you have skills to package software for this OS you are more than welcome to contribute by building those packages and / or writing documentation for these installs.

### III. Login

**Default Access :**

Login : admin

Password : admin



### IV. Orthanc Server declaration

By default the Orthanc server connexion settings is set to <http://localhost:8042>, with no login or password.  
If you have a different Orthanc location or if you add login / password, go to Options -> General and fill the form accordingly, click “update” to apply new connexion settings.  
Press “test connexion” to check that your Orthanc server can be correctly reached


# V. Main features

## a. Orthanc Content

The Orthanc Content tab, allows you to search for resources in your Orthanc Server.


You can use the search form to search for specific queries (done at Study level).

The results appears in dynamic table, the left column displays the Patient / Study level, at a study click the right column displays the series level.

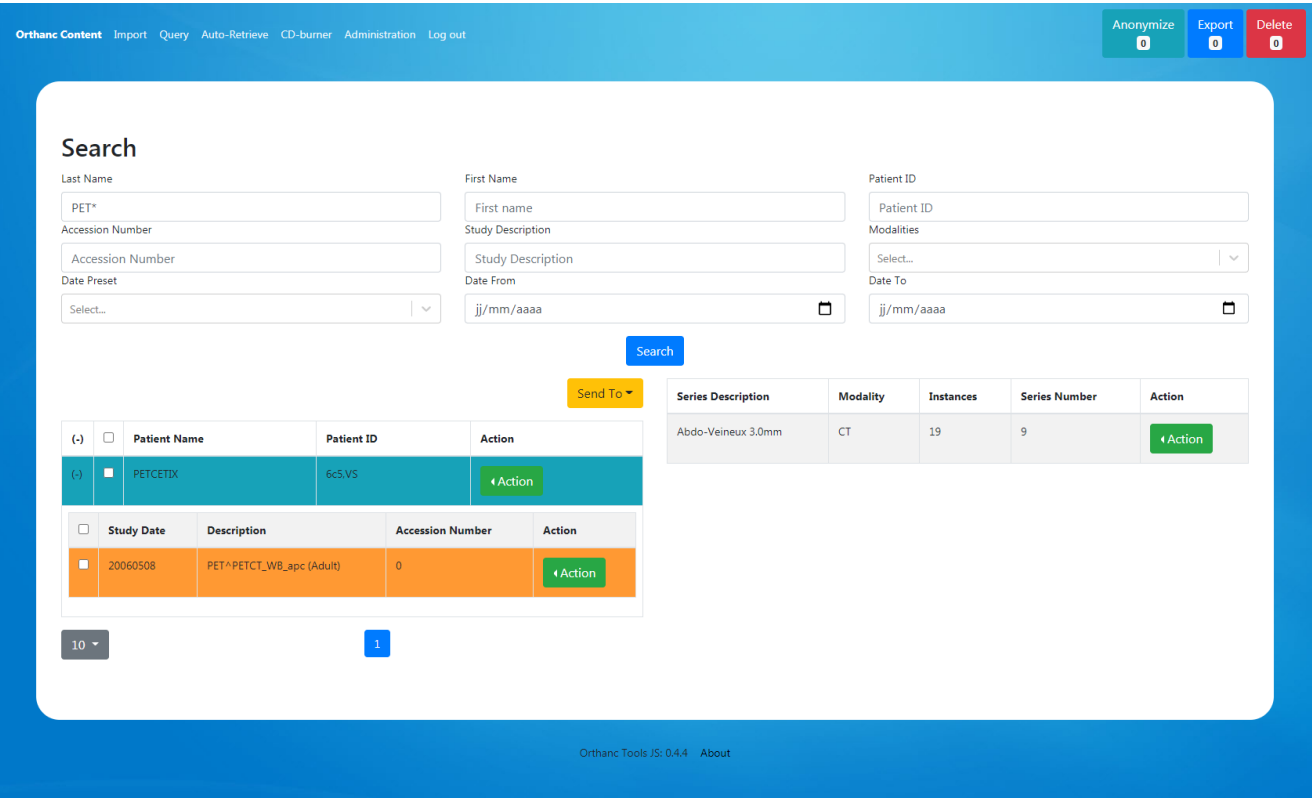

At each Level an “Action” button gives you access to additional functionalities such as:


- Modify (Dicom Edition, all levels)
- Delete (all levels)
- View on OHIF (Study level)
- Metadata (Series level)

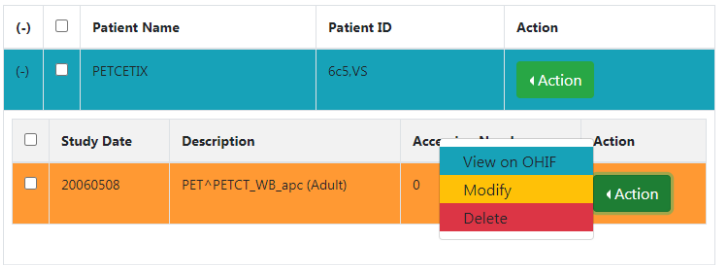

## ***b. Import Dicom files***

To import local dicom in Orthanc Server, simply drag and drop the dicoms files and folder in the dropzone.

The import is recursive in all folders and non dicoms files will not break import.

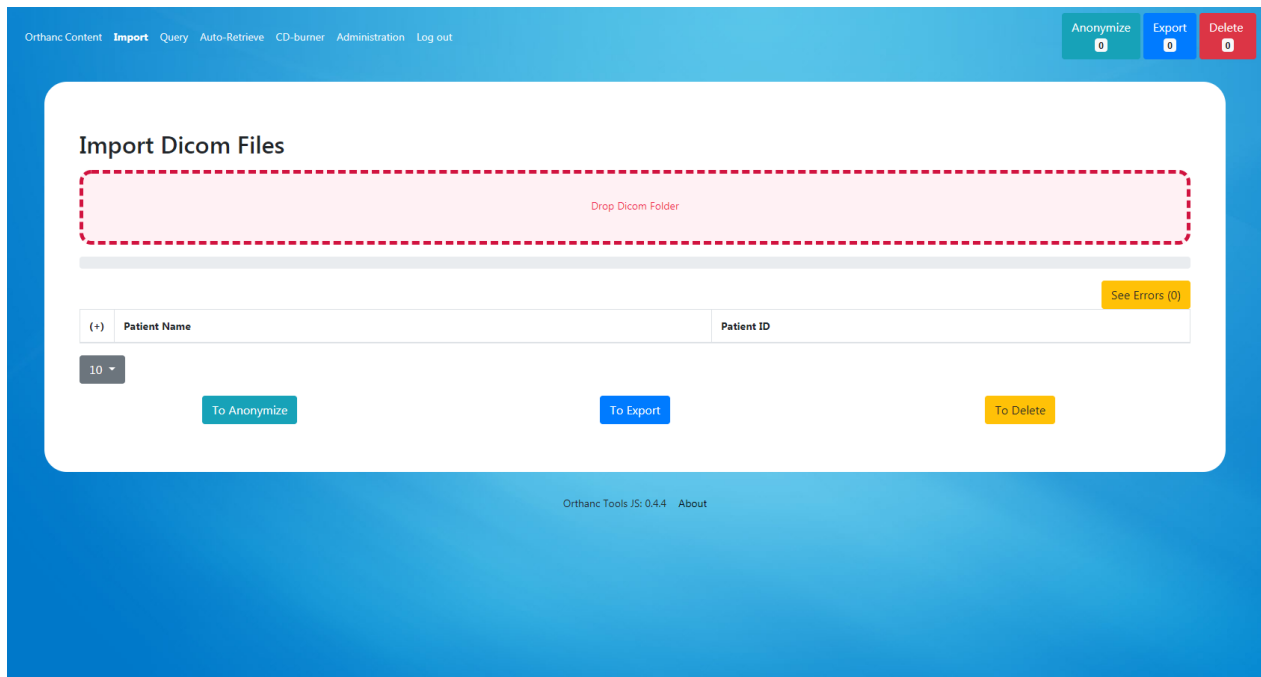

## ***c. Anonymize / Export / Delete***

At the right upper corner of Orthanc Tools JS you will see three tools: Anonymize, Export, Delete.

These tools can be seen as “baskets” to perform the anonymization, export or deletion action on a batch of resources.

You can add ressources (patients, studies) from many other modules of Orthanc Tools JS using the “To Anon List”, “To Export List”, “To Delete List” that you will find in several place of the software, allowing the dicom ressources to jump from a tool to another.

## ***i. Anonymize***

You can define a list of anonymisation, the anonymization is done at study level.

The GUI shows the patient level at the left table and study level in the right table.

The tables are editable to define, new patient name, patient ID, study description.

The Anonymization process follows the NEMA 2017c Guideline.

Two anonymization profiles are predefined in Orthanc Tools JS:

- Full Anonymization will follow the full Nema 2017c guideline
- Default Anonymization will keep body characteristics (weight, height...) and dates (acquisition date ...) this profile is recommended to anonymize nuclear medicine procedures (NM and PT).

Orthanc Content Import Query Auto-Retrieve CD-burner Administration Log out

Anonymize 1 Export 0 Delete 0

### Anonymize

Original Name	Original ID	New Name	New ID	Remove
PETCETIX	6c5.VS	<input type="text"/>	<input type="text"/>	<button>Remove</button>

10 1

prefix

Anon Profile :

Study Date	Description	New Description	Accession Number	New Accession Number	Remove
20060508	PET^PETCT_WB_apc (Adult)	PET^PETCT_WB_apc (Adult)	0	OrthancToolsJS	<button>Remove</button>

10 1

Orthanc Tools JS: 0.4.4 About

Once done, a result table will be shown allowing pushing anonymized resources to Export / Delete tool.

## ***ii. Export***

The export is a list of resources that can be export to:

- A simple browser download with transcoding capabilities
- To an AET
- To an Orthanc Peer
- To a FTP/SFTP/WebDav server (not yet ready)

Orthanc Content
Import
Query
Auto-Retrieve
CD-burner
Administration
Log out

Anonymize
0
Export
3
Delete
0

## Export

Patient Name	Patient ID	Study Date	Description	Accession Number	Anonymized ?
gfhgh0	gfhgh0	20180713		OrthancToolsJS	Yes
gfhgh0	gfhgh0	20180713		OrthancToolsJS	Yes
Femoral trombenarterectomy^Case Report:	Case Report 1	20110824			No

Series Description	Modality	Instances	Series Number	Remove
	MR	1200	11	Remove

10

1

Empty List

Download

Send To Modality

Send To Peer

Send To FTP

Send To WebDav

Implicit VR Endian

Orthanc Tools JS: 0.4.4 About

### iii. Delete

The deletion tool is a simple list a resource that will be deleted when the deletion process will be triggered.

Delete List

Empty List

(+)

Patient Name

Patient ID

Remove

(+)

Anonyme

3682

Remove

10

1

Delete List

Delete 1



## d. Query / Retrieve

Using a similar way to Orthanc Content you can do Query / Retrieve in your declared AETs.

The Study results are expandable to display related series.

Orthanc Content Import **Query** Auto-Retrieve CD-burner Administration Log out

Anonymize 0 Export 0 Delete 0

### Query

Last Name: g\* First Name: Patient ID: Accession Number: Study Description: Modalities: Select... Date Preset: Select... Date From: jj/mm/aaaa Date To: jj/mm/aaaa

PACS

Patient Name	Patient ID	Accession Number	Acquisition Date	Study Description	Modalities	Series	Instances	Retrieve
(-) Enter Patient Name	(-) Enter Patient ID...	(-) Enter Accession Nu	jj/mm/	(-) Enter Study Descri	(-) Enter Modalities...			
(-) GLIROPAD1	GLIROPAD1		20180713		MR	1	1200	Retrieve

Serie Description	Modality	Serie Number	Instances	Retrieve
	MR	11	1200	Retrieve

10 1

Orthanc Tools JS: 0.4.4 About

## e. Auto Retrieve


This feature allows you to make mass retrieve of DICOM from a PACS.

You can predefine a list of query to be done and retrieved and schedule it to be started outside working hour.


A robot will retrieve resources one by one to your Orthanc Server.

In the Query tab, define the list of queries you want to do:

- Click the "Add" button to open a new line and edit the wanted filed
- The table can be saved / imported using CSV file

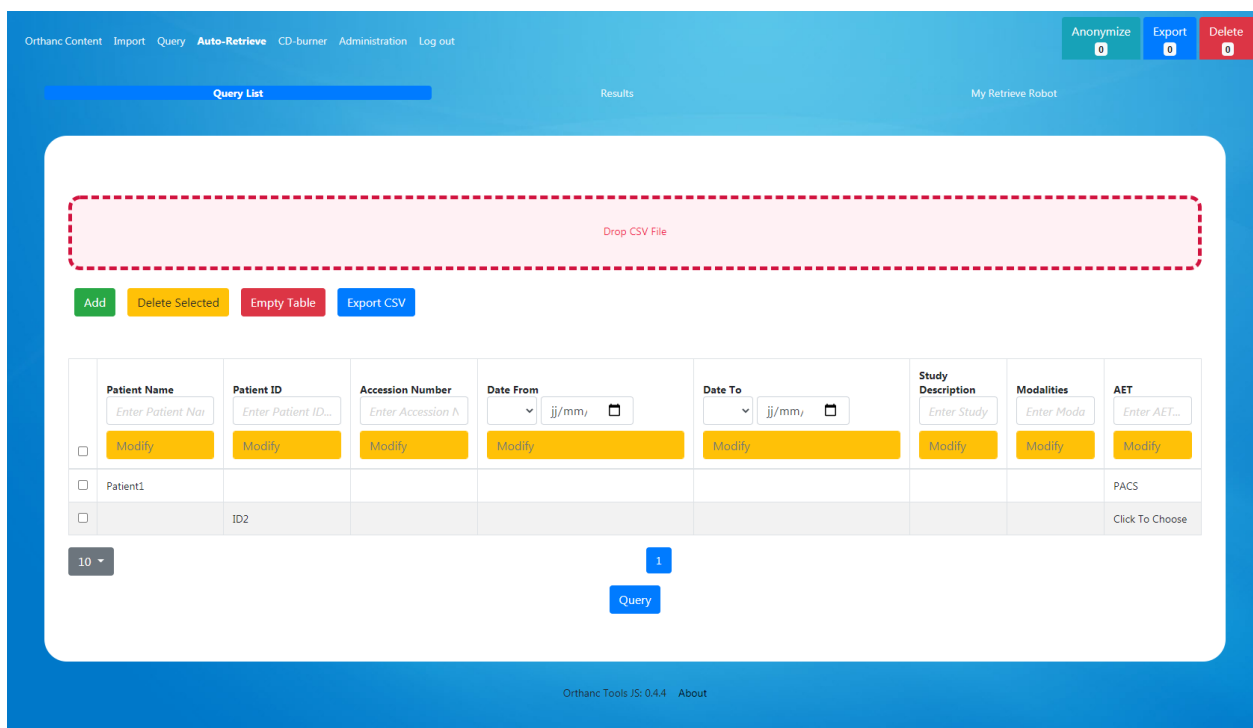

By clicking the “Query” Button, the defined queries will be done and results shown in the “results” tab.

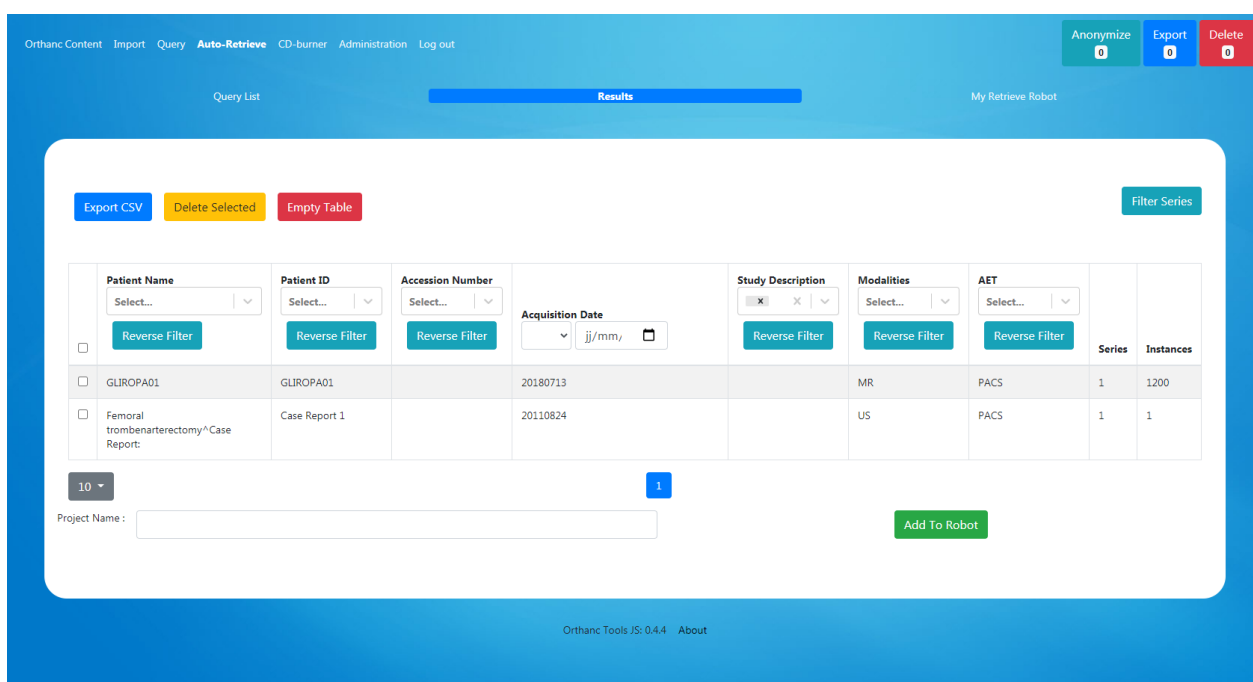

If you want to see series level click “filter series”, each query will be reprocessed at series level.

You can use the filters, search function of this table to search for your wanted resources. Remove all unwanted results, using the “delete selected” button.

At the bottom of this table you can define a project name and click “Add to Robot” the resources listed in the results table will be scheduled for automatic retrieve.

Orthanc Content
Import
Query
**Auto-Retrieve**
CD-burner
Administration
Log out

Anonymize
0
Export
0
Delete
0

Query List
Results
My Retrieve Robot

## Robot for user salim, project : myRobot

Progress : 0%

	level	Patient Name	Patient ID	Study Date	Modality	Study Description	Series Description	Accession Number	AET	Validated	Status	Remove Query	View in OHIF
		Enter Patient Name...	Enter F...	<div> <div></div> <div>jj/mm/yy</div> <div></div> </div>	Enter	Enter S...	Enter S...	Enter	E...	Enter	Ente...		
<input type="checkbox"/>	study	GLIOPA01	GLIOPA01	20180713	MR				PACS	false	Pending	Remove	View on OHIF
<input type="checkbox"/>	study	Femoral trombenarterectomy^Case Report	Case Report 1	20110824	US				PACS	false	Pending	Remove	View on OHIF

10
1

To Anonymize
To Export
To Delete

Orthanc Tools JS: 0.4.4 About

Your retrieve robot is scheduled. It needs to be validated by the administration (see administration panel).

It will start at the time set by the administrator.

Once retrieved the resources status will be marked "Retrieved" and you will be allowed to send retrieved resources to the anonymization / export tools of Orthanc Tools JS

## ***f. CD Burner – Epson / Primera DiscProducers***

This feature allows you to automatically trigger CD/DVD burning in commercial DiscProducers such as Epson (PP100, PP100N...) or Primera (SE, SE-3, II, Pro...)



The software can send burning request at the study (1 study per CD/DVD) or Patient level (multiple Study of a patient per CD/DVD). The burning request is automatically triggered when the dicom is stored in Orthanc (60 seconds after the last Dicom reception, see “StableAge” in Orthanc configuration if needed).

Along with Dicom files you can specify a viewer folder distribution to be associated with DICOM (such as Radiant for instance), and you can specify your own label to be printed in the CD/DVD which can include tags such as PatientName etc...

## ***i. Setting up burning solution***

To setup your burning solution, first go into the options tab and define:

- The Monitored Folder (Folder monitored by Epson or Primera application)
- The Viewer Folder containing viewer files to be included with Dicoms
- The Label File (see Epson/Primera template editor)
- The Transfer Syntax of Dicoms (if you want to burn a specific Dicom transfert syntax)
- The CD Burner Manufacturer (Epson or Primera)
- The monitoring Level (Study or Patient level)
- The Support Type (CD/DVD/Auto)
- If you want to automatically delete resources from Orthanc once sent to CD/Burner (to avoid infinity usage of disk space in Orthanc Server)

Orthanc Content Import Query Auto-Retrieve CD-burner **Administration** Log out

Anonymize 0 Export 0 Delete 0

General  
Aets  
Peers  
Robots  
Jobs  
CD Burner  
Plugins  
User Management

### CD/DVD Burner Options

Monitored Folder :  
C:\Users\kanoun\_s\Documents\cdBurner

Viewer Folder :  
C:\Users\kanoun\_s\Documents\Viewer

Label Path :  
C:\Users\kanoun\_s\Documents\label.ext

Transfer Syntax :  
JPEG Lossless

Manufacturer :  
Epson

Monitoring Level :  
Patient

Support Type :  
Auto

Delete Original Images From Orthanc ☒

Send

Orthanc Tools JS: 0.4.1 About

### ***For CD/DVD labels templates:***

Use Epson / Primera software to generate your CD/DVD printing and save it (.tdd file for Epson, .std for Primera) you can include the following tags (for Epson) or their respective position (for Primera)

“{patientName}” - Position 1

“{patientId}” – Position 2

“{studyDate}” - Position 3

“{studyDescription}” – Position 4

“ {patientDOB}” – Position 5

“ {accessionNumber}” – Position 6

to automatically print their respective value in the label.

### ***ii. Starting burning service***

Go do “CD-Burner” main tab and Start CD Burner Service, each time Orthanc will receive a Study / Patient, a burning request will be triggered and will appear in the interface (the burning service will continue to run even if the browser is closed)

Orthanc Content Import Query Auto-Retrieve **CD-burner** Administration Log out

Anonymize 0 Export 0 Delete 0

### CD Burner Service

☒ ☐

Patient Name	Patient ID	Patient Birth Date	Study Date	Study Description	CD Status	Cancel
10						

Orthanc Tools JS: 0.4.1 About

## VI. Administration of Orthanc-Tools-JS


### a. User Access

Two kinds of user definitions can be used


- Local users defined in the local database
- Distant users defined in an Active Directory

For each user a Role needs to be defined.

You can create how many Roles you want in the Role panel, each role is attached with a set of permission for each functionality of Orthanc-Tools-JS

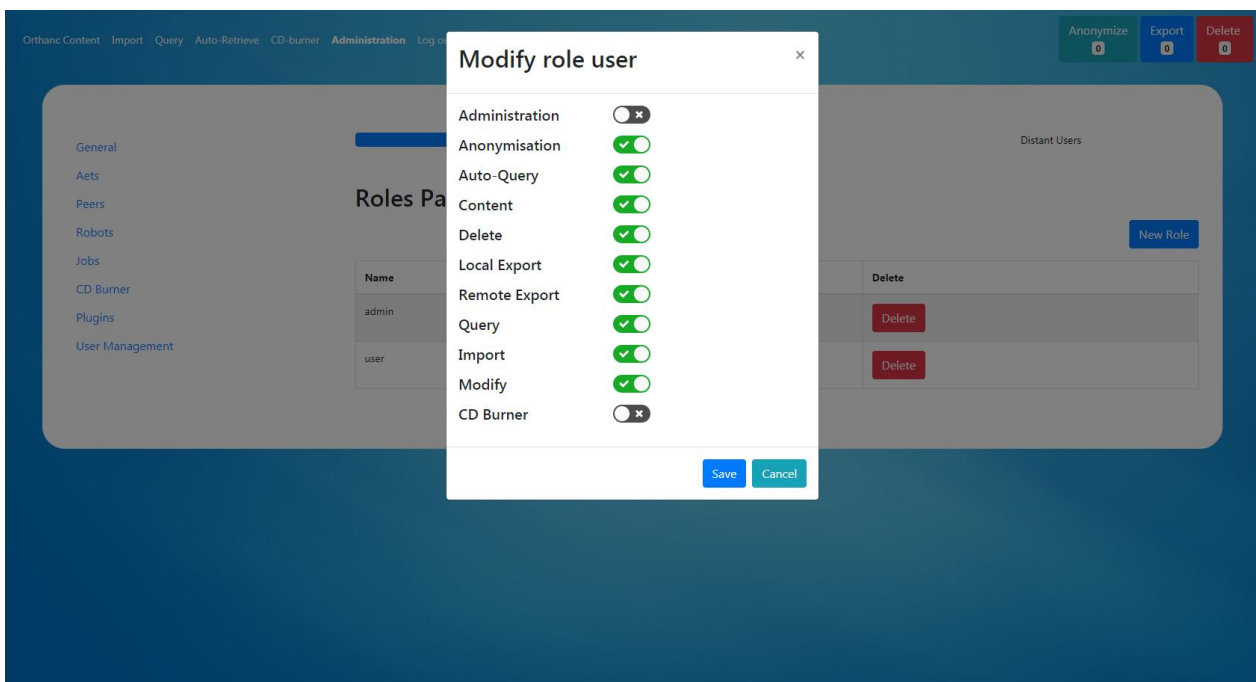


In the “Local User” tab you can define your local user, declared in database

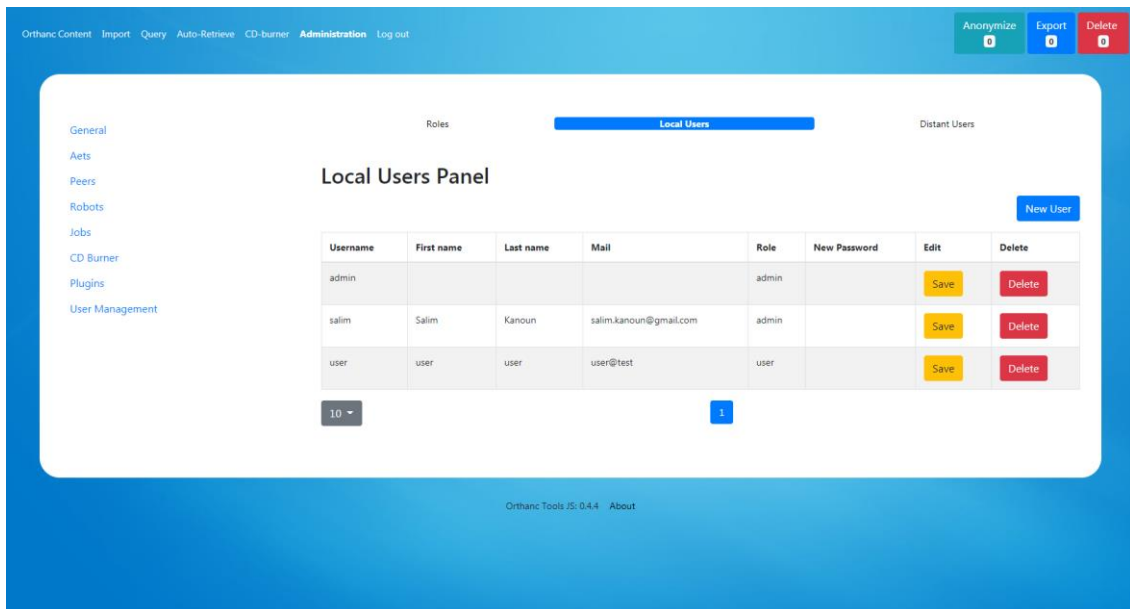


Last you can declare an Active Directory server, from with user will be granted access.

The Role definition is made by mapping the “Group Name” of the AD to one of your defined Role in Orthanc Tools JS

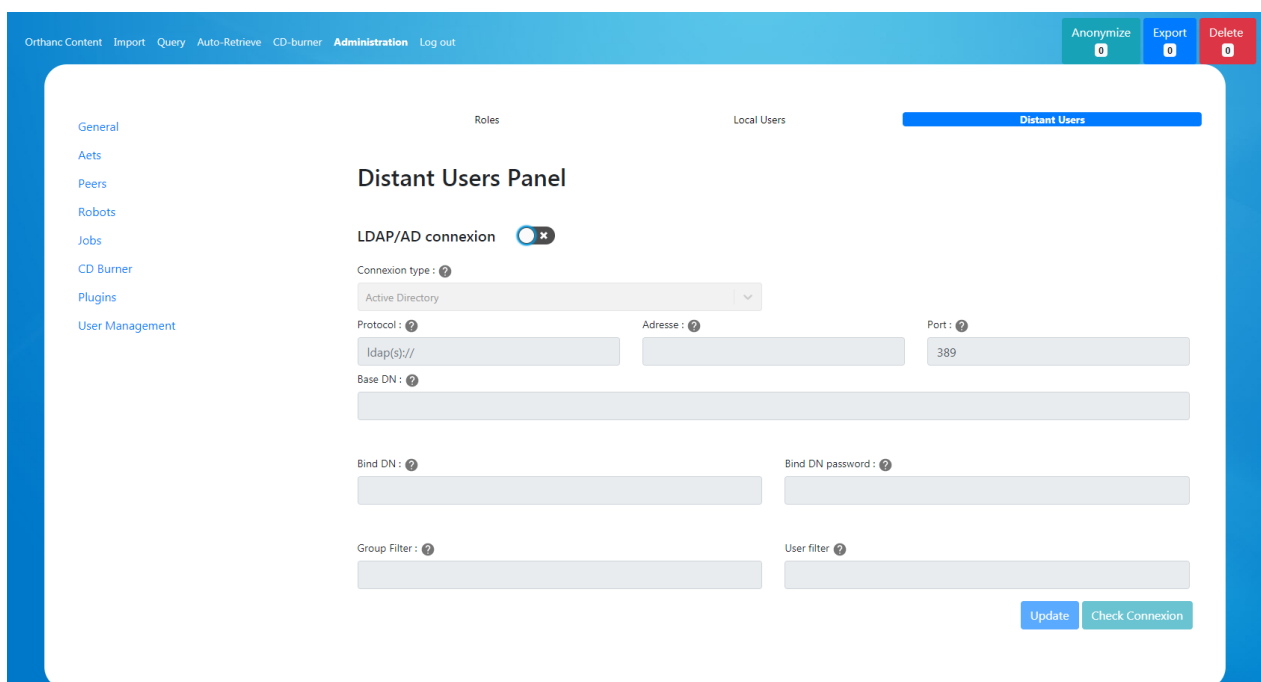

## ***b. AET / Orthanc Peers Declaration***

You can simply define / remove AET and Orthanc Peer in the Orthanc Server using the administration panel

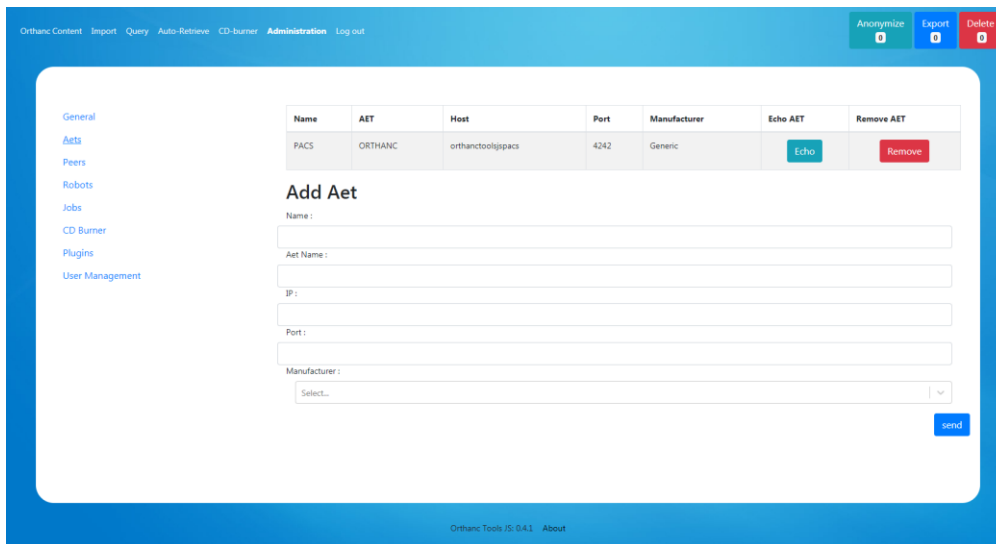

### c. Orthanc Jobs / Plugins

This interface allows you to see Orthanc asynchronous Jobs (and pause / cancel / resubmit them)

The plugin page list the available plugin in the current orthanc installation

### d. Auto Retrieve Robot

You can see robot request for automatic retrieve, asked by users.

The robot needs to be validated to be executed, during the validation, OrthancToolsJS will check that each query request provide only 1 answer in the PACS (to avoid too general queries that could lead to infinite retrieve).

You can see the robot query details and modify it if needed. You can also erase a retrieve job.

Last you can define the schedule hour in which the retrieve procedure will start (the retrieve are sequential)

